

AND

Digital Blood Pressure Monitor

Model UA-767NFC

Instruction Manual
Original

ENGLISH

Manuel d'instructions
Traduction

FRANÇAIS

Manual de Instrucciones
Traducción

ESPAÑOL

Manuale di Istruzioni
Traduzione

ITALIANO

Bedienungsanleitung
Übersetzung

DEUTSCH

Contents

| | |
|--|----|
| Dear Customers | 1 |
| Preliminary Remarks..... | 2 |
| Precautions..... | 2 |
| Parts Identification..... | 3 |
| Symbols | 4 |
| Operation Mode | 6 |
| Using the Monitor..... | 7 |
| Installing / Changing the Batteries..... | 7 |
| Connecting the AC Adapter | 7 |
| Connecting the Air Hose | 8 |
| Setting and Adjusting the Clock..... | 8 |
| Notes on the Built-in Clock..... | 9 |
| Selecting the Correct Cuff Size..... | 10 |
| Applying the Arm Cuff | 11 |
| How to Take Proper Measurements | 11 |
| Measurement | 11 |
| After Measurement | 11 |
| Measurements..... | 12 |
| Normal Measurement..... | 12 |
| Measurement with the Desired Systolic Pressure..... | 13 |
| Near Field Communication (NFC) Function | 15 |
| Displaying the Memory Data..... | 15 |
| Switching the Memory Setting..... | 16 |
| Transmitting the Memory Data | 17 |
| Transmitting the Measurement Data | 17 |
| What is an Irregular Heartbeat..... | 18 |
| Pressure Bar Indicator..... | 18 |
| About Blood Pressure | 18 |
| What is Blood Pressure?..... | 18 |
| What is Hypertension and How is it Controlled? | 19 |
| Why Measure Blood Pressure at Home? | 19 |
| WHO Blood Pressure Classification | 19 |
| Blood Pressure Variations..... | 19 |
| Troubleshooting..... | 20 |
| Maintenance | 21 |
| Technical Data..... | 21 |

Dear Customers

Congratulations on purchasing a state-of-the-art A&D blood pressure monitor, one of the most advanced monitors available today. Designed for ease of use and accuracy, this device will facilitate your daily blood pressure regimen.

We recommend that you read through this manual carefully before using the device for the first time.

Preliminary Remarks

- This device conforms to the European Directive 93/42 EEC for Medical Products. This is made evident by the **CE**₀₁₂₃ mark of conformity.
(0123: The reference number to the involved notified body)
- The device is designed for use on adults only, not newborns or infants.
- Environment for use. The device is for use in the home healthcare environment.
- This device is designed to measure blood pressure and pulse rate of people for diagnosis.

Precautions

- Precision components are used in this device. Extreme temperature, humidity, direct sunlight, shock or dust should be avoided.
- Clean the device and cuff with a dry, soft cloth or a cloth dampened with water and a neutral detergent. Never use alcohol, benzene, thinner or other harsh chemicals to clean the device or cuff.
- Avoid tightly folding the cuff or storing the hose tightly twisted for long periods, as such treatment may shorten the life of the components.
- The device and cuff are not water resistant. Prevent rain, sweat and water from soiling the device and cuff.
- Measurements may be affected if the device is used near a television, microwave oven, cellular telephone, X-ray or other devices that emit electromagnetic waves.
- Used equipment, parts and batteries are not treated as ordinary household waste, and must be disposed of according to the applicable local regulations.
- When the AC adapter is used, make sure that the AC adapter can be readily removed from the electrical outlet when necessary.
- When reusing the device, confirm that the device is clean.
- Do not modify the device. It may cause accidents or damage to the device.
- To measure blood pressure, the arm must be squeezed by the cuff hard enough to temporarily stop blood flow through the artery. This may cause pain, numbness or a temporary red mark to the arm. This condition will appear especially when measurement is repeated successively. Any pain, numbness, or red marks will disappear with time.


Contraindications

The following are precautions for proper use of the device.


- Do not apply the cuff on an arm with another medical electrical equipment attached. The equipment may not function properly.
- People who have a severe circulatory deficit in the arm must consult a doctor before using the device, to avoid medical problems.
- Do not self-diagnose the measurement results and start treatment by yourself. Always consult your doctor for evaluation of the results and treatment.
- Do not apply the cuff on an arm with an unhealed wound.
- Do not apply the cuff on an arm receiving an intravenous drip or blood transfusion. It may cause injury or accidents.

- Do not use the device where flammable gases such as anesthetic gases are present. It may cause an explosion.
- Do not use the device in highly concentrated oxygen environments, such as a high-pressure oxygen chamber or an oxygen tent. It may cause a fire or explosion.

Parts Identification

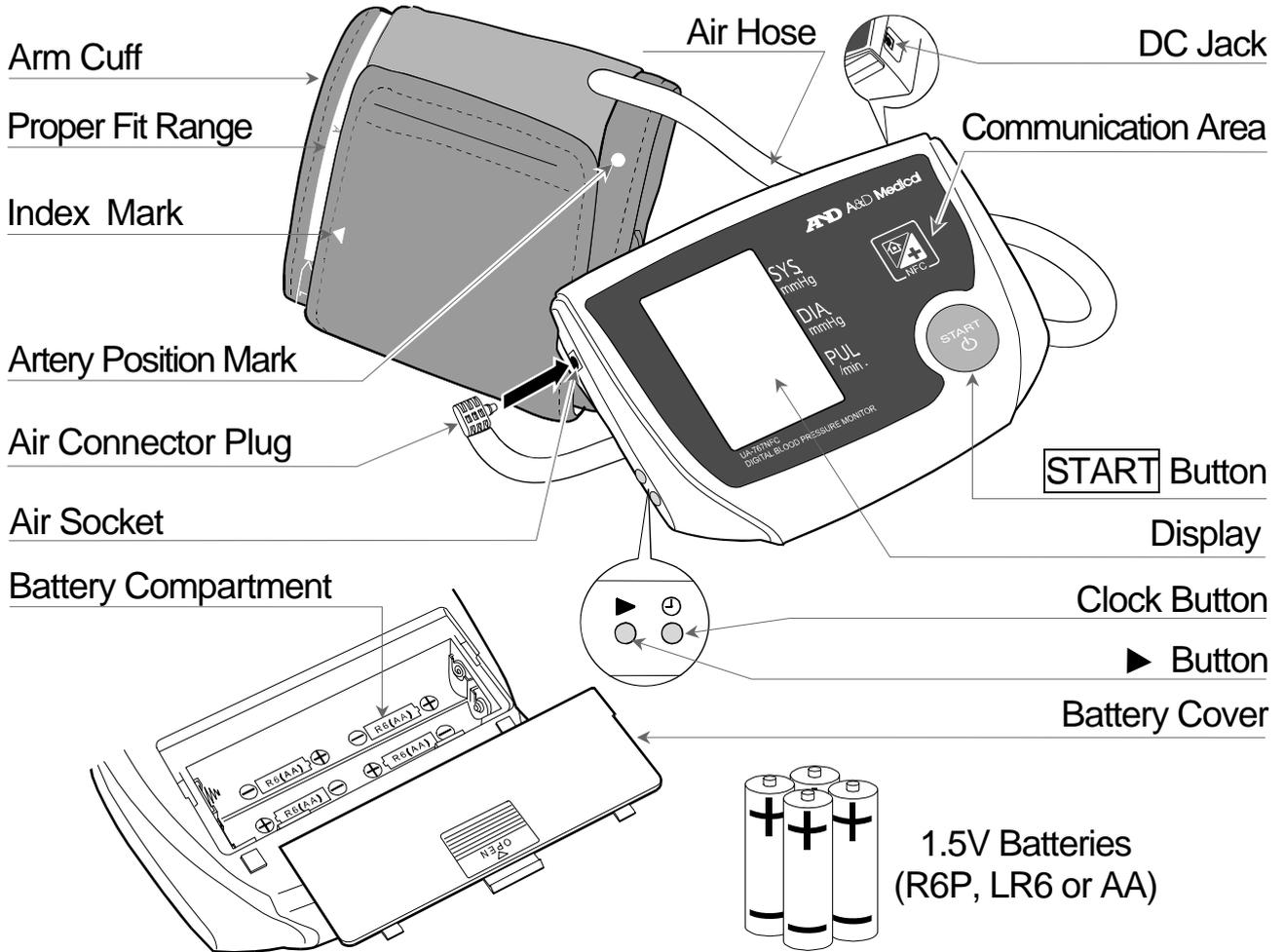


Part of Display

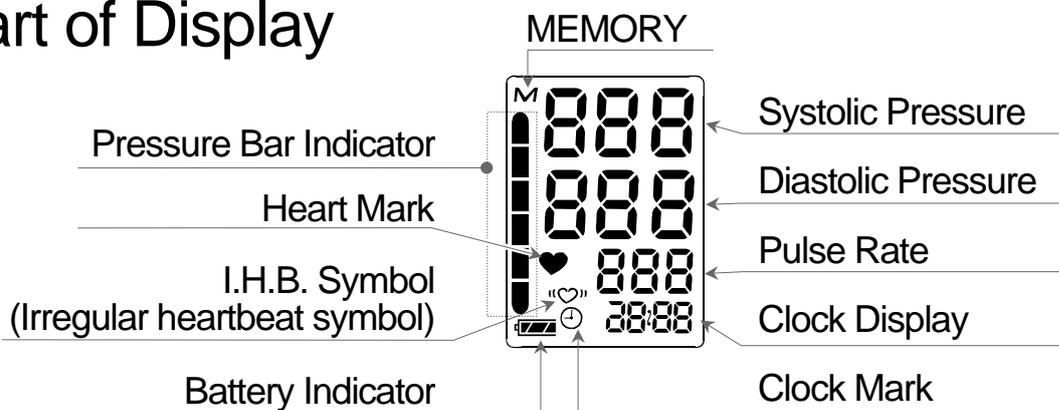

Symbols

Symbols that are printed on the device case
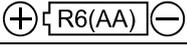
| Symbols | Function / Meaning |
|---|--|
|  | Standby and Turn the device on. |
| SYS. | Systolic blood pressure in mmHg |
| DIA. | Diastolic blood pressure in mmHg |
| PUL./min. | Pulse per minute |
|  | Battery installation guide |
|  | Direct current |
|  | Type BF: Device, cuff and tubing are designed to provide special protection against electrical shocks. |
|  | EC directive medical device label |
|  | EU-representative |
|  | Manufacturer |
|  | Date of manufacture |
|  | Class II device |
|  | WEEE label |
| SN | Serial number |
|  | Refer to instruction manual/booklet |
|  | Polarity of DC jack |

Symbols that appear on the display

| Symbols | Function / Meaning | Recommended Action |
|---|---|--|
|  | Appears while measurement is in progress. It blinks when the pulse is detected. | Measurement is in progress. Remain as still as possible. |
|  | Irregular Heartbeat symbol (I.H.B.) Appears when an irregular heartbeat is detected. It may light when a very slight vibration like shivering or shaking is detected. | _____ |
| M | Previous measurements stored in memory. | _____ |

Symbols

| Symbols | Function / Meaning | Recommended Action |
|---|---|--|
|  Full Battery | The battery power indicator during measurement. | _____ |
|  Low Battery | The battery is low when it blinks. | Replace all batteries with new ones when the mark blinks. |
| Err | Unstable blood pressure due to movement during measurement. | Take another measurement. Remain very still during measurement. |
| | The systolic and diastolic values are within 10 mmHg of each other. | Apply the cuff correctly, and take another measurement. |
| | The pressure value did not increase during the inflation. | |
|  Err CUF | The cuff is not applied correctly. | |
|  Err PUL DISPLAY ERROR | The pulse is not detected correctly. | |
|  Err 10 | Cannot communicate. | The communication device may be slanted or placed out of the communication area of the device. Hold the communication symbol of the NFC-enabled communication device close to the communication area of the device until "End" appears. |
|  Err 11 | | |
|  Err 12 | | |
|  Err 13 | Attempted to transmit the memory data when the memory setting is 0. | Set the memory to 100. |
|  | Appears during the clock setting procedure | _____ |

Operation Mode

1. Normal Measurement

Press the **START** button. Blood pressure is measured.

2. Memory Function

When the memory setting is 100:

This device will store the last 100 measurements in memory.
When the measurement exceeds 100, the oldest data will be deleted automatically and the most recent data will be stored.

When the memory setting is 0:

This device will not store any data in memory.

3. Deleting all Data Stored in Memory


Make sure that the device displays the clock.

Press and hold the **▶** button.

The **M** mark appears. After a while it starts blinking, then disappears.

Now all data stored in memory has been deleted.

Note: This operation will delete all data stored in memory.
You cannot select which data to delete.

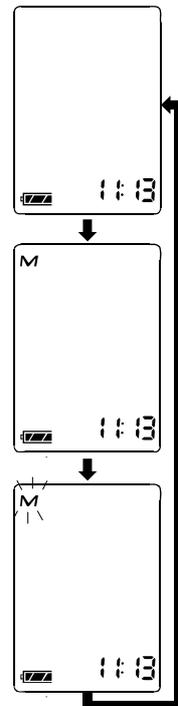


4. Measurement with the Desired Systolic Pressure

Refer to page 13 for measurement with the desired systolic pressure.

Using the Monitor

Installing / Changing the Batteries

1. Slide the battery cover up to open it.
2. Remove the used batteries and insert new batteries into the battery compartment as shown, taking care that the polarities (+ and -) are correct. Use only R6P, LR6 or AA batteries.
3. Slide the battery cover down to close.

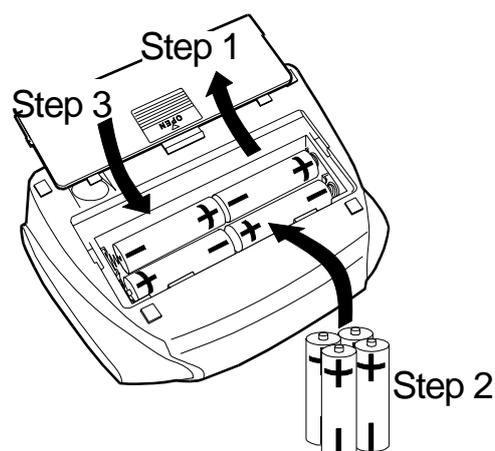


Connecting the AC Adapter

Insert the AC adapter plug into the DC jack.

Next, connect the AC adapter to an electrical outlet.

- Use the specified AC adapter. (Refer to page 22.)

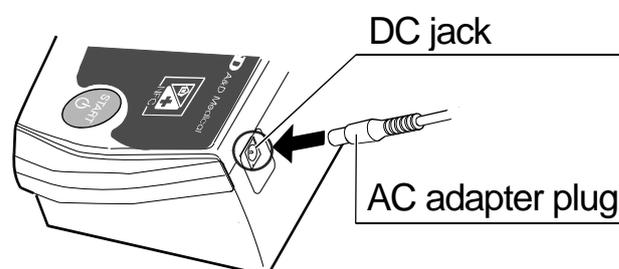

Note: When inserting the batteries or connecting the AC adapter, “Year” blinks automatically and the screen will go to the clock setting mode. Set the clock as described on page 8.

Every time the batteries are inserted or the AC adapter is connected, the clock needs to be reset. Measurement cannot start until a time is set on the clock.


CAUTION

- Insert the batteries as shown in the battery compartment. If installed incorrectly, the device will not work.
- When  (LOW BATTERY mark) blinks on the display, replace all batteries with new ones. Do not mix old and new batteries. It may shorten the battery life, or cause the device to malfunction. Wait two seconds or more after turning the device off, to replace the batteries. If  (LOW BATTERY mark) appears even after the batteries are replaced, make a blood pressure measurement. The device may then recognize the new batteries.
-  (LOW BATTERY mark) does not appear when the batteries are drained.
- The battery life varies with the ambient temperature and may be shorter at low temperatures.
- Use the specified batteries only.
- Remove the batteries if the device is not to be used for a long time. The batteries may leak and cause a malfunction.

Using the Monitor

Connecting the Air Hose

Insert the air connector plug into the air socket firmly.

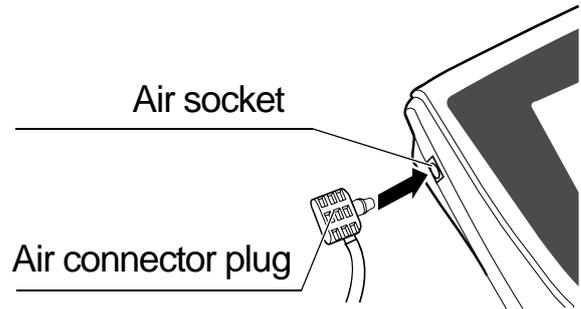

Setting and Adjusting the Clock


Set or adjust the clock prior to use.

1. Press and hold the  button until the year starts blinking.

When inserting the batteries or connecting the AC adapter, "Year" blinks automatically.

2. Select the year using the  button. Press the  button to set the current year and move to month/day selection. The date can be set anywhere between the years 2011 and 2055.
3. Select the month using the  button. Press the  button to set the current month and move to day selection.
4. Select the day using the  button. Press the  button to set the current day and move to hour/minute selection.
5. Select the hour using the  button. Press the  button to set the current hour and move to minute selection.
6. Select the minute using the  button. Press the  button to set the current minute and activate the clock.

Press and hold  "Year" blinks when inserting the batteries or connecting the AC adapter

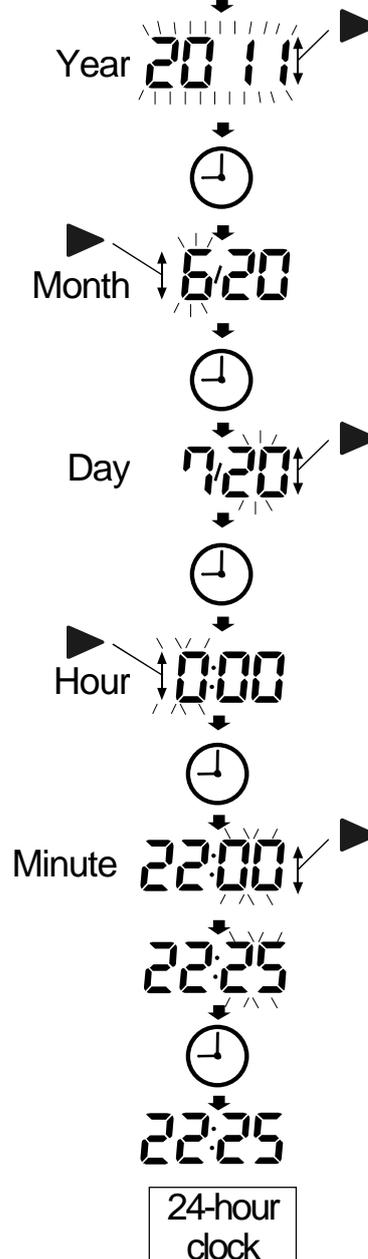

Notes on the Built-in Clock

- The device displays the clock as long as power is connected.
- Adjust the clock prior to use.
- During measurement, the clock adjustment is not available.
- When adjusting the clock, press the ► button to increase the value by one, or press and hold the ► button to change the value continuously.
- The setting procedure cannot be reversed. If the wrong value is set, start the setting procedure again.
- The current clock setting will be cleared if the batteries are removed, or power is not connected for approximately 30 seconds while the AC adapter is used. In this case, adjust the clock again.

If the clock is set and you want to manually adjust the time:

- Press and hold the ⌚ button to go to clock settings.
- You can cancel the settings by pressing the [START] button and the screen will return to the clock display mode.
- After 30 seconds of non-operation, the setting procedure will be automatically canceled and the screen will return to the clock display mode.

When inserting, exchanging the batteries or connecting the AC adapter:

- “Year” blinks automatically and the screen will go to clock settings.
- To start measurement, set the clock as described on page 8.
- Pressing the [START] button will not return to the clock display mode. You must complete clock settings first.
- After 30 seconds of non-operation, the setting procedure will be automatically canceled and the screen will display 0:00. (Clock setting is not completed and measurement cannot start.)
- When exchanging the batteries or reconnecting the AC adapter, the screen displays Year, month, and date last time the batteries or AC adapter was removed. You need to complete the setting of actual date and clock before use.

Using the Monitor

Selecting the Correct Cuff Size

Using the correct cuff size is important for an accurate reading. If the cuff is not the proper size, the reading may yield an incorrect blood pressure value.

- The arm size is printed on each cuff.
- The index ▲ and proper fit range, on the cuff, tell you if you are applying the correct cuff. (Refer to "Symbols that are printed on the cuff" below)
- If the index ▲ points outside of the range, contact your local dealer to purchase a replacement cuff.
- The arm cuff is a consumable. If it becomes worn, purchase a new one.


| Arm Size | Recommended Cuff Size | Catalog Number |
|----------------|-----------------------|----------------|
| 31 cm to 45 cm | Large adult cuff | CUF-F-LA |
| 22 cm to 32 cm | Adult Cuff | CUF-F-A |

Arm size: The circumference of the biceps.

Note: Model UA-767NFC is NOT designed for using a small adult cuff.

Symbols that are printed on the cuff.


| Symbols | Function / Meaning | Recommended Action |
|------------|--|---|
| ● | Artery Position Mark | Set the ● mark on the artery of the upper arm or in line with the ring finger on the inside of the arm. |
| ▲ | Index | _____ |
| REF | Catalog number | _____ |
| A | Proper fit range for the Adult cuff. It's printed on the Adult cuff. | _____ |
| L | Range to use the Large adult cuff. Over range printed on the Adult cuff. | Use the Large adult cuff instead of the Adult cuff. |
| L | Proper fit range for the Large adult cuff. It's printed on the Large adult cuff. | _____ |
| S | Under range printed on the Adult cuff. | |
| A | Range to use the Adult cuff. It's printed on the Large adult cuff. | Use the Adult cuff instead of the Large adult cuff. |
| LOT | Lot number | _____ |

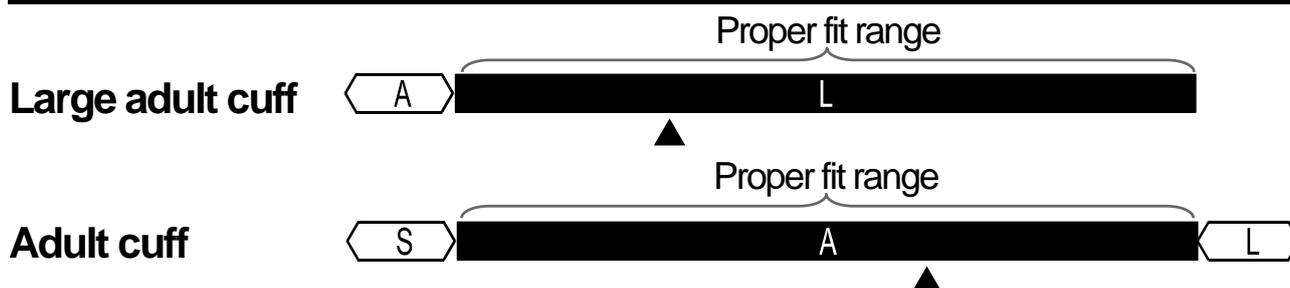

Using the Monitor

Applying the Arm Cuff

1. Wrap the cuff around the upper arm, about 1 – 2 cm above the inside of the elbow, as shown. Place the cuff directly against the skin, as clothing may cause a faint pulse, and result in a measurement error.
2. Constriction of the upper arm, caused by tightly rolling up a shirtsleeve, may prevent accurate readings.
3. Confirm that the index ▲ points within the proper fit range.

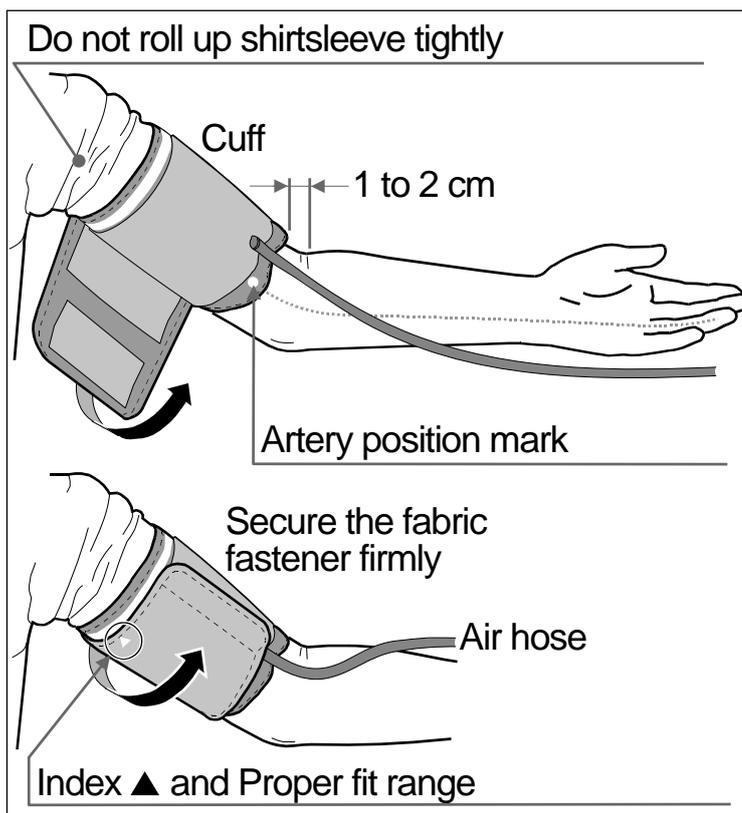

How to Take Proper Measurements

For the most accurate blood pressure measurement:

- Sit comfortably on a chair. Rest your arm on the table. Do not cross your legs. Keep your feet on the floor and straighten your back.
- Relax for about five to ten minutes before measurement.
- Place the center of the cuff at the same height as your heart.
- Remain still and keep quiet during measurement.
- Do not measure right after physical exercise or a bath. Rest for twenty or thirty minutes before taking the measurement.
- Try to measure your blood pressure at the same time every day.

Measurement

During measurement, it is normal for the cuff to feel very tight. (Do not be alarmed).

After Measurement

After measurement, press the **START** button to turn the device off. After one minute of non-operation, the device will turn off automatically. Remove the cuff and record your data.

Note: Allow at least three minutes between measurements on the same person.

Measurements

Before measurement, read “Notes for Proper Measurement” on page 14.

Normal Measurement

1. Place the cuff on the arm (preferably the left arm).
Sit quietly during measurement.

2. Press the **START** button.
All of the display segments are displayed.
Zero (0) is displayed blinking briefly.

The display changes, as indicated in the figure at the right, as the measurement begins. The cuff starts to inflate. It is normal for the cuff to feel very tight. A pressure bar indicator is displayed, on the left edge of the display, during the inflation.

Note: If you wish to stop inflation at any time, press the **START** button again.

3. When inflation is complete, deflation starts automatically and ♥ (heart mark) blinks, indicating that the measurement is in progress. Once the pulse is detected, the mark blinks with each pulse beat.


Note: If an appropriate pressure is not obtained, the device starts to inflate again automatically.

To avoid re-inflation, see “Measurement with the Desired Systolic Pressure” on the next page.

4. When the measurement is complete, the systolic and diastolic pressure readings and pulse rate are displayed.
The cuff exhausts the remaining air and deflates completely.

5. Press the **START** button to turn the device off.
After one minute of non-operation, the device will turn off automatically.

Note: Allow at least three minutes between measurements on the same person.

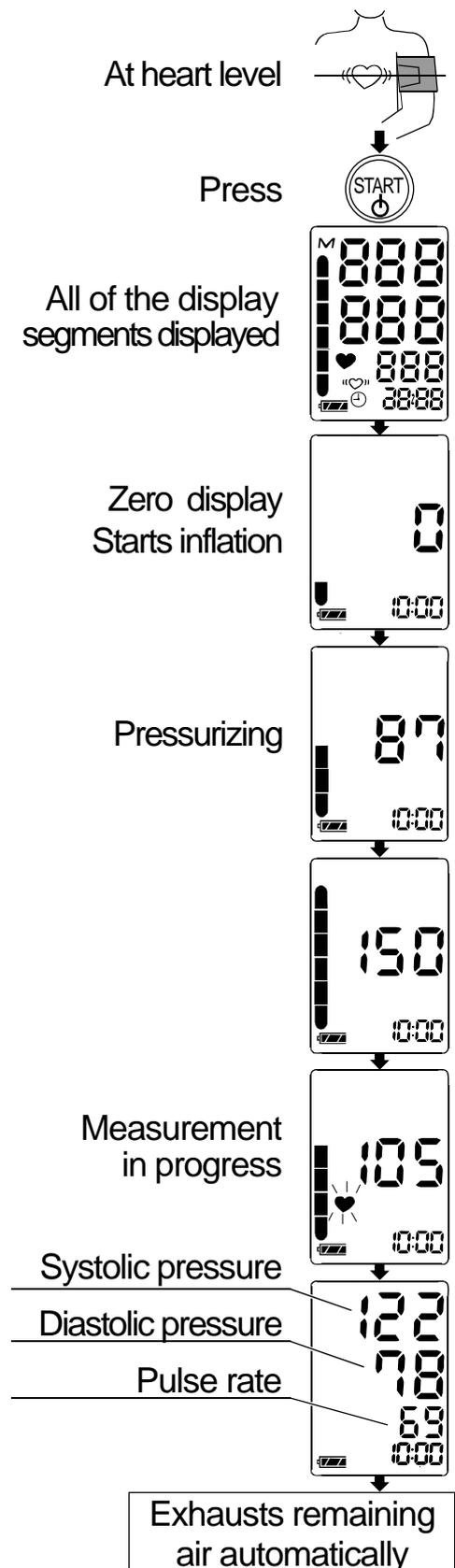

Measurements

Measurement with the Desired Systolic Pressure


Model UA-767NFC is designed to detect the pulse and to inflate the cuff to a systolic pressure level automatically.

Use this method when re-inflation occurs repeatedly or when the results are not displayed even if the pressure decreases to 20 mmHg or less.

When the memory setting is 0:

1. Place the cuff on the arm (preferably the left arm).
2. Press and hold the **START** button until a number about 30 to 40 mmHg higher than your expected systolic pressure appears.
3. When the desired number is reached, release the **START** button to start measurement. Continue to measure your blood pressure as described on the previous page.


At heart level

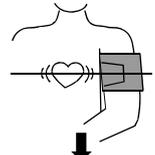

Press and hold the button



Release the button at the desired systolic pressure

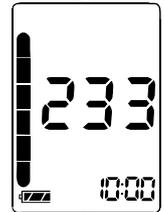


Refer to the previous page for measurement

When the memory setting is 100:

1. Place the cuff on the arm (preferably the left arm).
2. Press the **START** button. All of the display segments are displayed. Zero (0) is displayed blinking briefly.
3. Press and hold the **START** button until a number about 30 to 40 mmHg higher than your expected systolic pressure appears.
4. When the desired number is reached, release the **START** button to start measurement. Continue to measure your blood pressure as described on the previous page.


At heart level

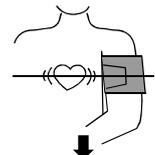


Press



All of the display segments displayed

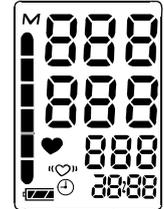

Zero blinking

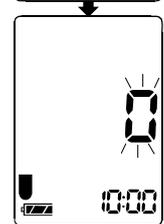

Press and hold the button



Release the button at the desired systolic pressure



Refer to the previous page for measurement

Notes for Proper Measurement

- Sit down in a comfortable position. Place your arm on a table with your palm facing upward and the cuff at the same level as your heart.
- Relax for about five to ten minutes before taking a measurement. If you are excited or depressed by emotional stress, the measurement will reflect this stress as a higher (or lower) than normal blood pressure reading and the pulse reading will usually be faster than normal.
- An individual's blood pressure varies constantly, depending on what you are doing and what you have eaten. What you drink can have a very strong and rapid effect on your blood pressure.
- This device bases its measurements on the heartbeat. If you have a very weak or irregular heartbeat, the device may have difficulty determining your blood pressure.
- Should the device detect a condition that is abnormal, it will stop the measurement and display an error symbol. Refer to page 5 for the description of symbols.
- This device is intended for use by adults only. Consult with your physician before using this device on a child. A child should not use this device unattended.

Near Field Communication (NFC) Function

To transmit data to or from the device, communication devices that comply with the communication specifications of the device are required.

Examples of communication devices that comply

- NFC-enabled cellular phone
- Personal computer with NFC reader/writer
- Access Point


Displaying the Memory Data

When the memory setting is 100, the memory data can be displayed.

While the device displays the clock, press and hold the **START** button. The systolic pressure, diastolic pressure and pulse rate with the memory number are displayed. At every 1.5 seconds, the memory number display changes to the date, then to the time, and then to the memory number again.

Press the  button to increase the memory number by one and press the  button to decrease the memory number by one.

Press the **START** button to return to the clock display. After one minute of non-operation, the device will return to the clock display.

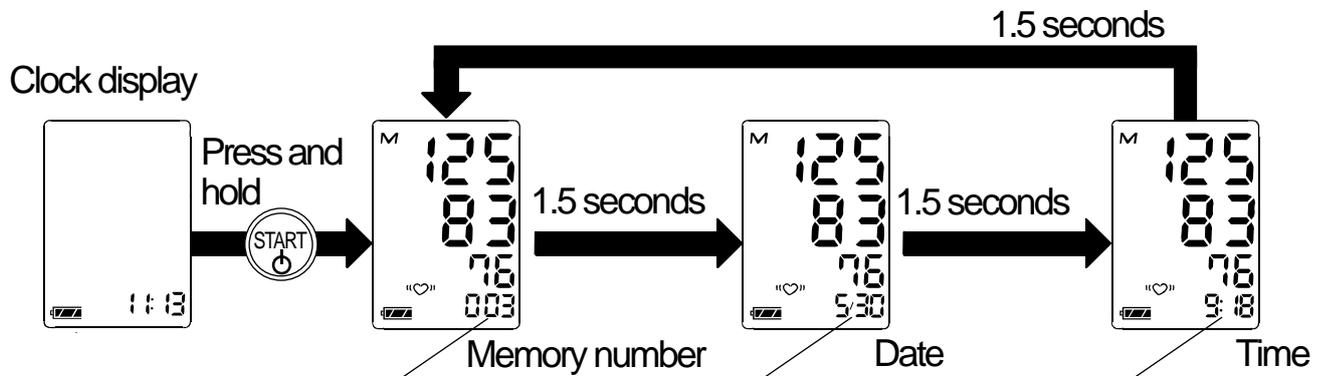

Press to increase the memory number by one.



Press to decrease the memory number by one.





Press to return to the clock display.

Note: Displaying the memory data is not available when the memory setting is 0.

Switching the Memory Setting


While the device displays the clock, hold the communication symbol of the NFC-enabled communication device close to the communication area () of the device until “End” appears. The memory setting is switched.

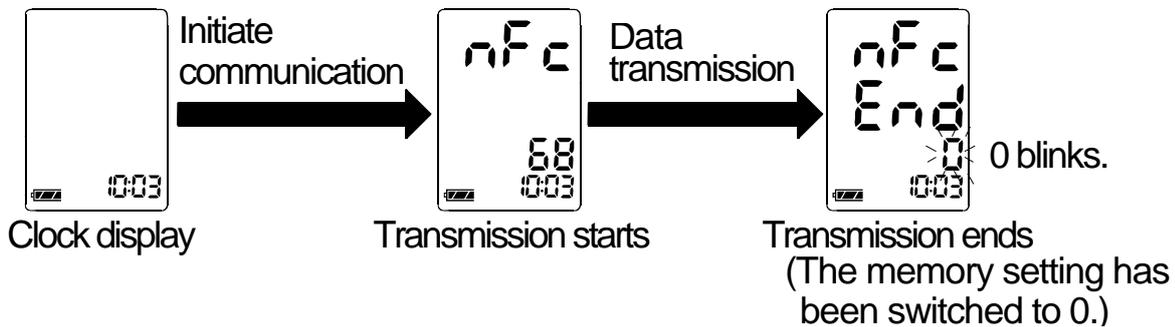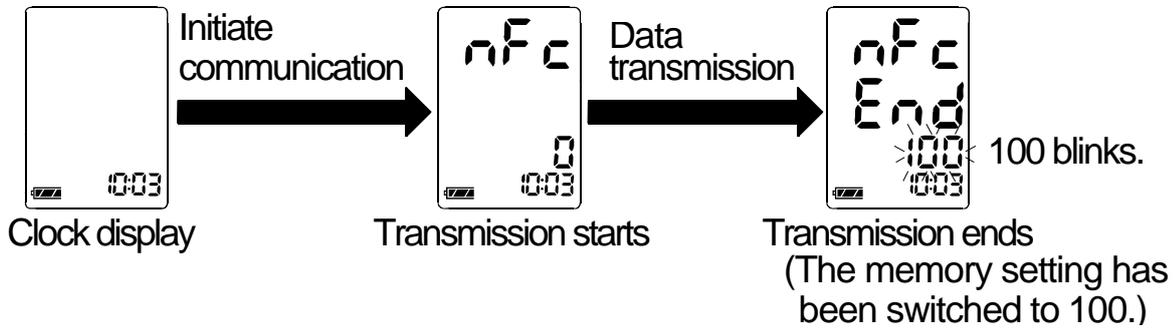

Note: When switching the memory setting from 100 to 0, the memory data, if any, will be deleted.


Transmitting the Memory Data

While the device displays the clock, hold the communication symbol of the NFC-enabled communication device close to the communication area (📶) of the device until “End” appears. The memory data transmission is performed. If the clock of the device is not correct, the clock will be updated at the same time.

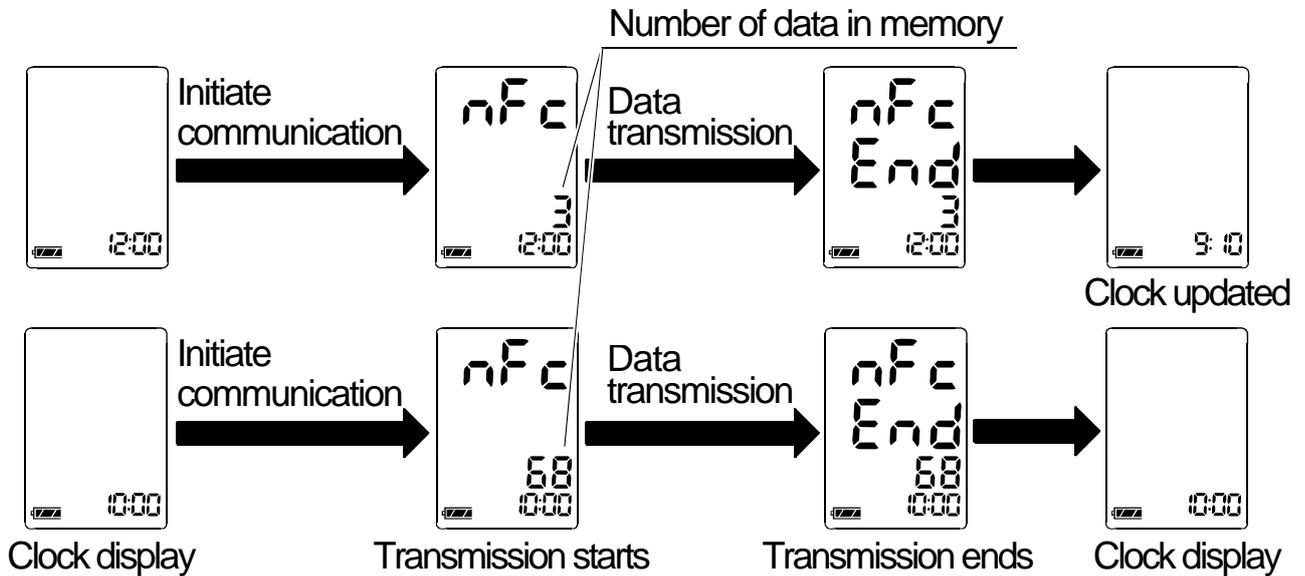

Transmitting the Measurement Data

When the measurement is complete and the device displays the results, hold the communication symbol of the NFC-enabled communication device close to the communication area (📶) of the device until “End” appears. The data transmission is performed. If the clock of the device is not correct, the clock will be updated at the same time.

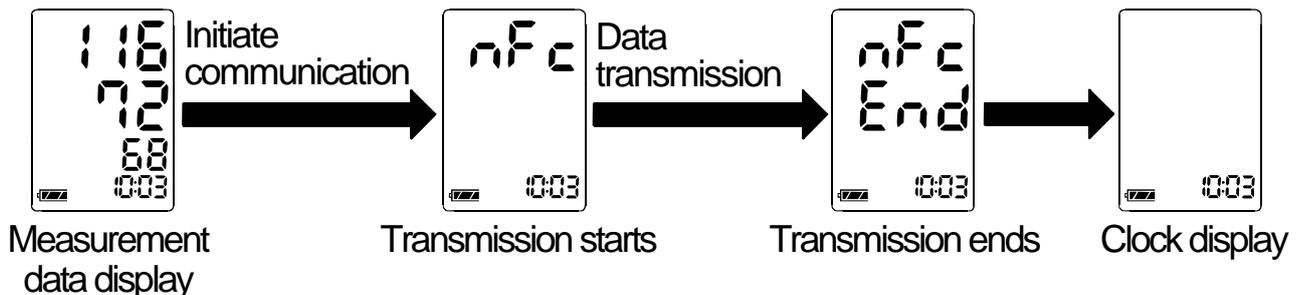

Note: The result display will turn off automatically after one minute. To transmit the measurement data, initiate the communication within one minute after the measurement is complete.

Caution on Using the NFC Function

The UA-767NFC has an NFC Forum Type 3 Tag wireless interface module built in. When the device is used near wireless communication devices which use the same frequency as that of the device (13.56 MHz) or use its harmonic frequencies, electromagnetic interference may occur. In that case, turn off the devices that are not in use, or place the device 1 meter or more away from the devices.


What is an Irregular Heartbeat

The UA-767NFC provides a blood pressure and pulse rate measurement even when an irregular heartbeat occurs. An irregular heartbeat is defined as a heartbeat that varies by 25% from the average of all heartbeats during the blood pressure measurement. It is important that you are relaxed, remain still and do not talk during measurements.

Note: We recommend contacting your physician if you see this (♥) indicator frequently.

Pressure Bar Indicator

The indicator monitors the progress of pressure during measurement.

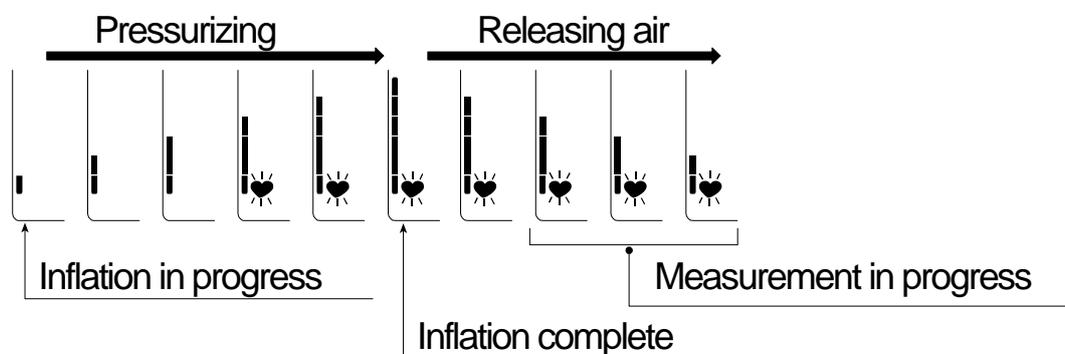

About Blood Pressure

What is Blood Pressure?

Blood pressure is the force exerted by blood against the walls of the arteries. Systolic pressure occurs when the heart contracts. Diastolic pressure occurs when the heart expands. Blood pressure is measured in millimeters of mercury (mmHg). One's natural blood pressure is represented by the fundamental pressure, which is measured first thing in the morning while one is still at rest and before eating.

About Blood Pressure

What is Hypertension and How is it Controlled?

Hypertension, an abnormally high arterial blood pressure, if left unattended, can cause many health problems including stroke and heart attack. Hypertension can be controlled by altering lifestyle, avoiding stress, and with medication under a doctor's supervision.

To prevent hypertension or keep it under control:

- Do not smoke
- Reduce salt and fat intake
- Maintain proper weight
- Exercise regularly
- Have regular physical checkups


Why Measure Blood Pressure at Home?

Blood pressure measured at a clinic or doctor's office may cause apprehension and can produce an elevated reading, 25 to 30 mmHg higher than that measured at home. Home measurement reduces the effects of outside influences on blood pressure readings, supplements the doctor's readings and provides a more accurate, complete blood pressure history.

WHO Blood Pressure Classification

Standards to assess high blood pressure, without regard to age, have been established by the World Health Organization (WHO), as shown in the chart.


Reference Material: Journal of Hypertension 1999, Vol 17 No.2

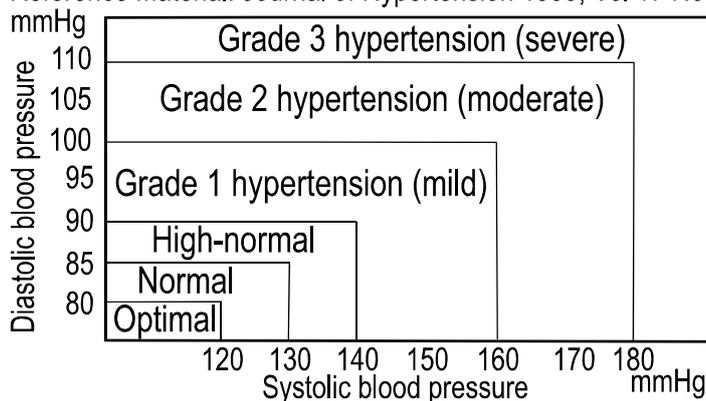

Blood Pressure Variations

An individual's blood pressure varies greatly on a daily and seasonal basis. It may vary by 30 to 50 mmHg due to various conditions during the day. In hypertensive individuals, variations are even more pronounced. Normally, the blood pressure rises while at work or play and falls to its lowest levels during sleep. So, do not be overly concerned by the results of one measurement.

Take measurements at the same time every day using the procedure described in this manual to get to know your normal blood pressure. Regular readings give a more comprehensive blood pressure history. Be sure to note the date and time when recording your blood pressure. Consult your doctor to interpret your blood pressure data.

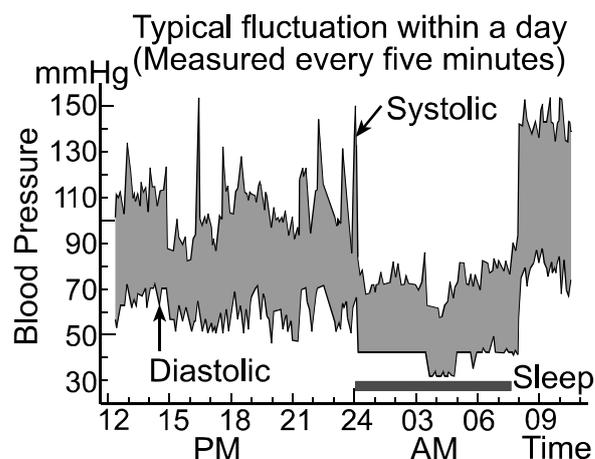

Troubleshooting

| Problem | Possible Reason | Recommended Action |
|---|---|---|
| Nothing appears on the display, even when the power is turned on. | Batteries are drained. | Replace all batteries with new ones. |
| | Battery terminals are not in the correct position. | Reinstall the batteries with negative and positive terminals matching those indicated on the battery compartment. |
| The cuff does not inflate. | Battery voltage is too low.  (LOW BATTERY mark) blinks. If the batteries are drained completely, the mark does not appear. | Replace all batteries with new ones. |
| The device does not measure. Readings are too high or too low. | The cuff is not applied properly. | Apply the cuff correctly. |
| | You moved your arm or body during measurement. | Make sure you remain very still and quiet during measurement. |
| | The cuff position is not correct. | Sit comfortably and still. Place your arm on a table with your palm facing upward and the cuff at the same level as your heart. |
| | _____ | If you have a very weak or irregular heart beat, the device may have difficulty in determining your blood pressure. |
| Other | The value is different from that measured at a clinic or doctor's office. | Refer to "Why Measure Blood Pressure at Home?". |
| | _____ | Remove the batteries. Place them back properly and take another measurement. |

Note: If the actions described above do not solve the problem, contact the dealer. Do not attempt to open or repair this product, as any attempt to do so will make your warranty invalid.

Maintenance

Do not open the device. It uses delicate electrical components and an intricate air unit that could be damaged. If you cannot fix the problem using the troubleshooting instructions, contact the authorized dealer in your area or our customer service department. The A&D customer service will provide technical information, spare parts and units to authorized dealers.

The device was designed and manufactured for a long service life. However it is generally recommended to have the device inspected every 2 years, to ensure proper functioning and accuracy. Please contact the authorized dealer in your area or A&D for maintenance.

Technical Data

| | |
|-------------------------------|--|
| Type | UA-767NFC |
| Measurement method | Oscillometric measurement |
| Measurement range | Pressure: 0 - 280 mmHg Pulse: 40 - 180 beats / minute |
| Measurement accuracy | Pressure: ± 3 mmHg Pulse: $\pm 5\%$ |
| Power supply | 4 x 1.5V batteries (R6P, LR6 or AA) or AC adapter (TB-233) (Not included) |
| Number of measurements | Approximately 450 measurements, when AA Alkaline batteries are used, with pressure value of 180 mmHg at room temperature of 23 °C. |
| Classification | Internally powered ME equipment (Supplied by batteries) / Class II (Supplied by adapter) Continuous operation mode |
| Clinical test | According to ANSI / AAMI SP-10 1992 |
| EMC | IEC 60601-1-2: 2007 |
| Memory | 0 or 100 measurements (Only for devices that rewriting of the memory data is enabled.) |
| Operating condition | +10 °C to +40 °C / 15 %RH to 85 %RH / 800 hPa to 1060 hPa |
| Transport / Storage condition | -10 °C to +60 °C / 15 %RH to 95 %RH |
| Dimensions | Approx. 147 [W] x 64 [H] x 110 [D] mm |
| Weight | Approx. 300 g, excluding the batteries |
| Applied part | Cuff Type BF  |

| | |
|------------------------|---|
| Useful life | Device: 5 years (when used six times a day) Cuff: 2 years (when used six times a day) |
| Wireless communication | NFC wireless interface module |
| Accessory AC adapter | The adapter is to connect the device to a power source at home. |
| TB-233 | Please contact your local A&D dealer for purchasing. The AC adapter is required to be inspected or replaced periodically. |
| TB-233C | Input: 100-240V Output: 6V  500mA      139°C 2A |
| TB-233BF | Input: 240V Output: 6V  500mA      139°C 2A |

Accessories sold separately

Cuff

| Catalog Number | Cuff Size | Arm Size |
|----------------|------------------|----------------|
| CUF-F-LA | Large adult cuff | 31 cm to 45 cm |
| CUF-F-A | Adult cuff | 22 cm to 32 cm |

Arm size: The circumference at the biceps.

AC adapter

| Catalog Number | Plug (Outlet type) |
|----------------|--------------------|
| TB-233C | Type C |
| TB-233BF | Type BF |

Note: Specifications are subject to change without prior notice.

EMC table information is listed on our website:

http://www.aandd.jp/products/manual/medical/emc_acadapter_en.pdf



 **A&D Company, Limited**

1-243 Asahi , Kitamoto-shi, Saitama 364-8585, JAPAN
Telephone: [81] (48) 593-1111 Fax: [81] (48) 593-1119

| | |
|----|-----|
| EC | REP |
|----|-----|

A&D INSTRUMENTS LIMITED

Unit 24/26 Blacklands Way, Abingdon Business Park, Abingdon, Oxfordshire OX14 1DY United Kingdom
Telephone: [44] (1235) 550420 Fax: [44] (1235) 550485

A&D ENGINEERING, INC.

1756 Automation Parkway, San Jose, California 95131, U.S.A.
Telephone: [1] (408) 263-5333 Fax: [1] (408)263-0119

A&D AUSTRALASIA PTY LTD

32 Dew Street, Thebarton, South Australia 5031, AUSTRALIA
Telephone: [61] (8) 8301-8100 Fax: [61] (8) 8352-7409